

Workstation Critical Care / Neonatal Care / EvitaV300 Option GS500




The GS500 serve as....:

- pneumatically backup in case of CGS failure
- continuous gas supply for hospitals without Central Gas Supply
- gas supply for inner clinical transport

The GS500 is capable to:

- drive the workplace for ca. 120min (with PS500)
- substitute the old bulky Compressor on workplaces with V500 / VN500

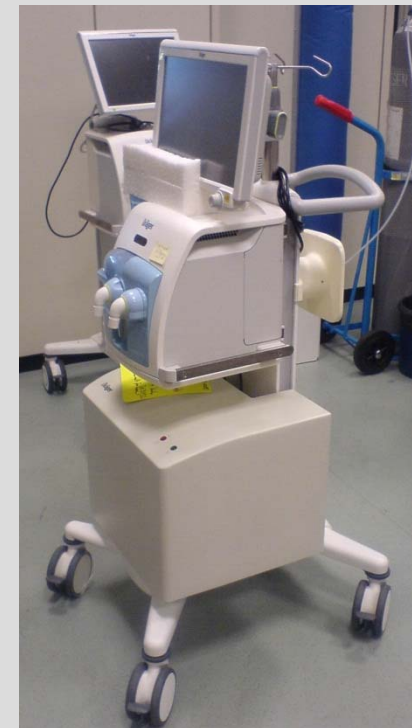

Workstation Critical Care / Neonatal Care / EvitaV300 Option GS500




- no user interaction necessary-
automatically switch on/off
- usable for continuous gas delivery or
as gas-back-up-system
- GS500 and PS500 support inner-
clinical transport.
- low noise
- fixed mounted – no additional cable
or hoses disturb user
- HFO is not possible during transport
with GS500!

Workstation Critical Care / Neonatal Care / EvitaV300 Option GS500




Gas delivery from GS500 to V-unit is running via the gas-channel, which is integrated into column of trolley.

- **Transport function**
 - start up GS500
 - via transport-button in dialog „special maneuver“
 - or
 - disconnect the hoses from central gas supply

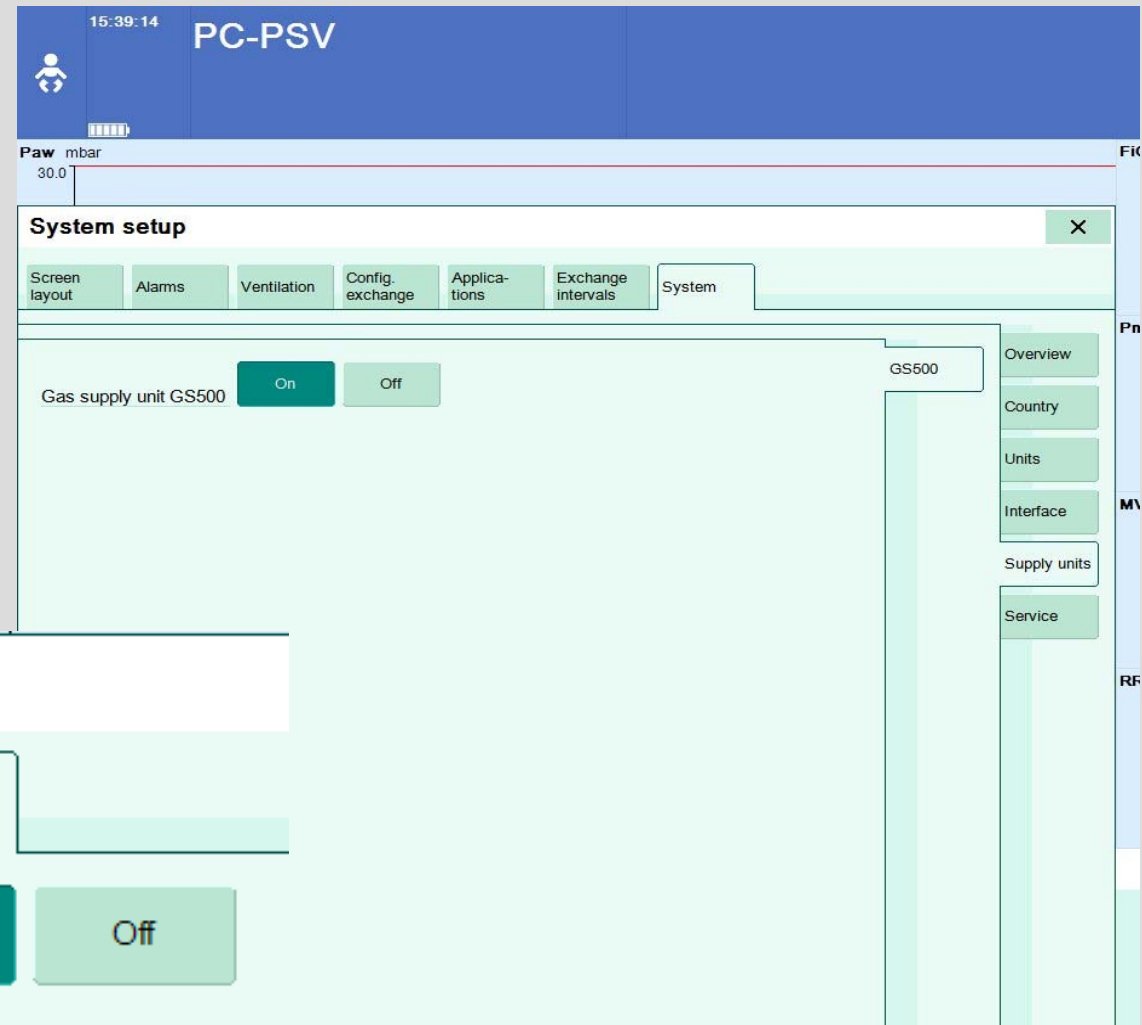

Workstation Critical Care / Neonatal Care / EvitaV300

Option GS500

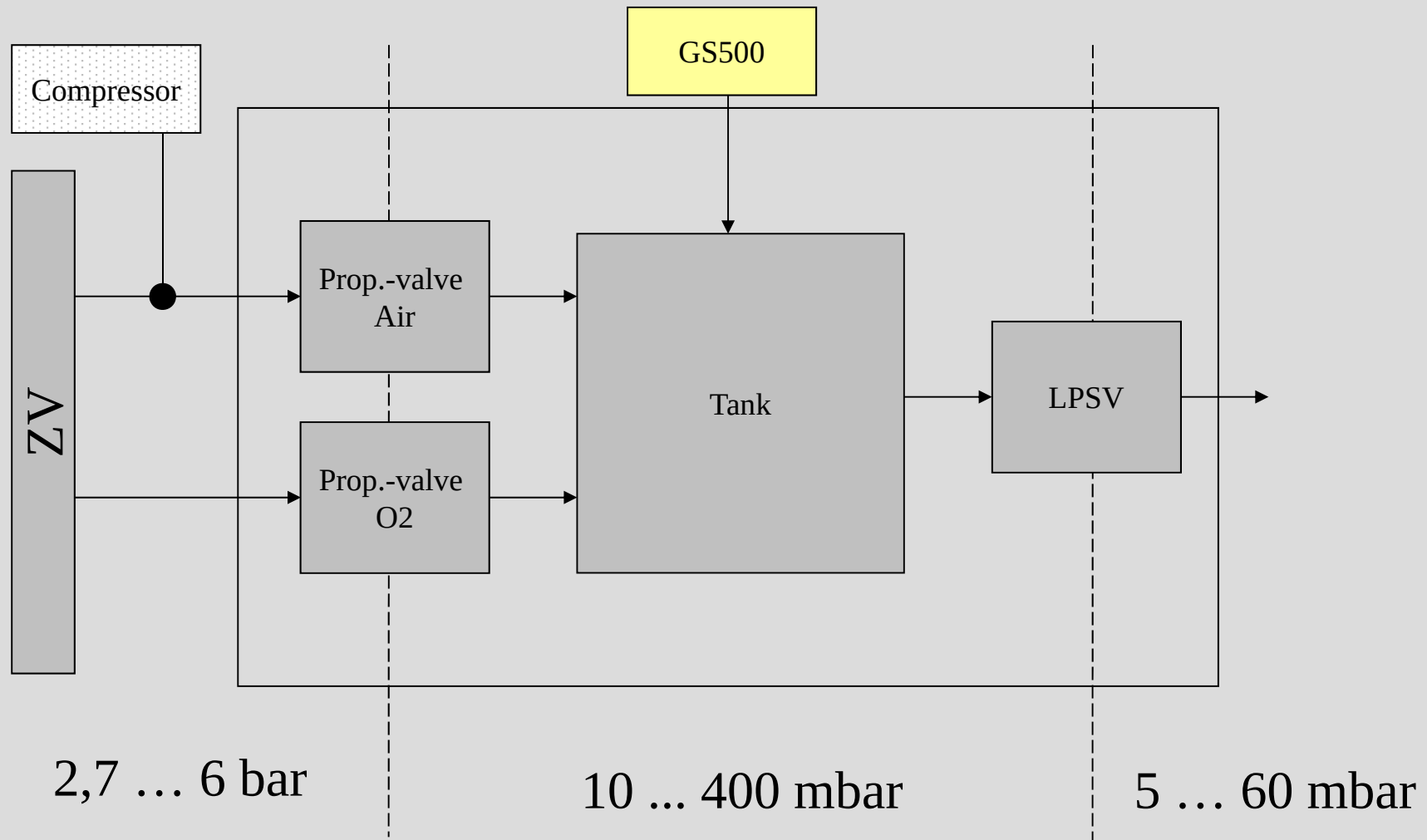


2009-12-01/ C500/V500 Training, G. Hagelstein

V-Unit with GS500


Tankpressure approx. 20mbar higher than P_{insp}

Inspiration

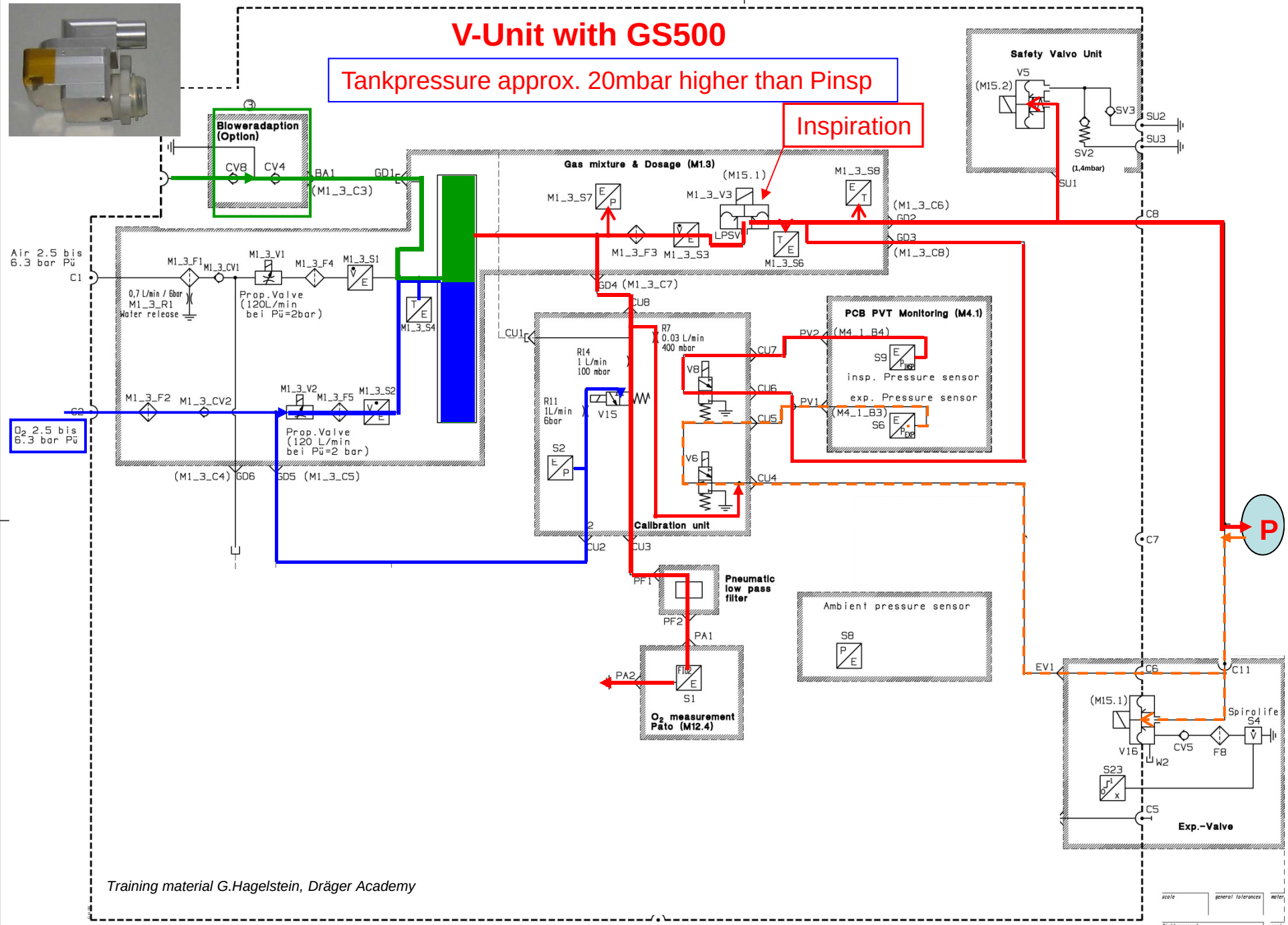


Workstation Critical Care / Neonatal Care / EvitaV300

Option GS500, function description

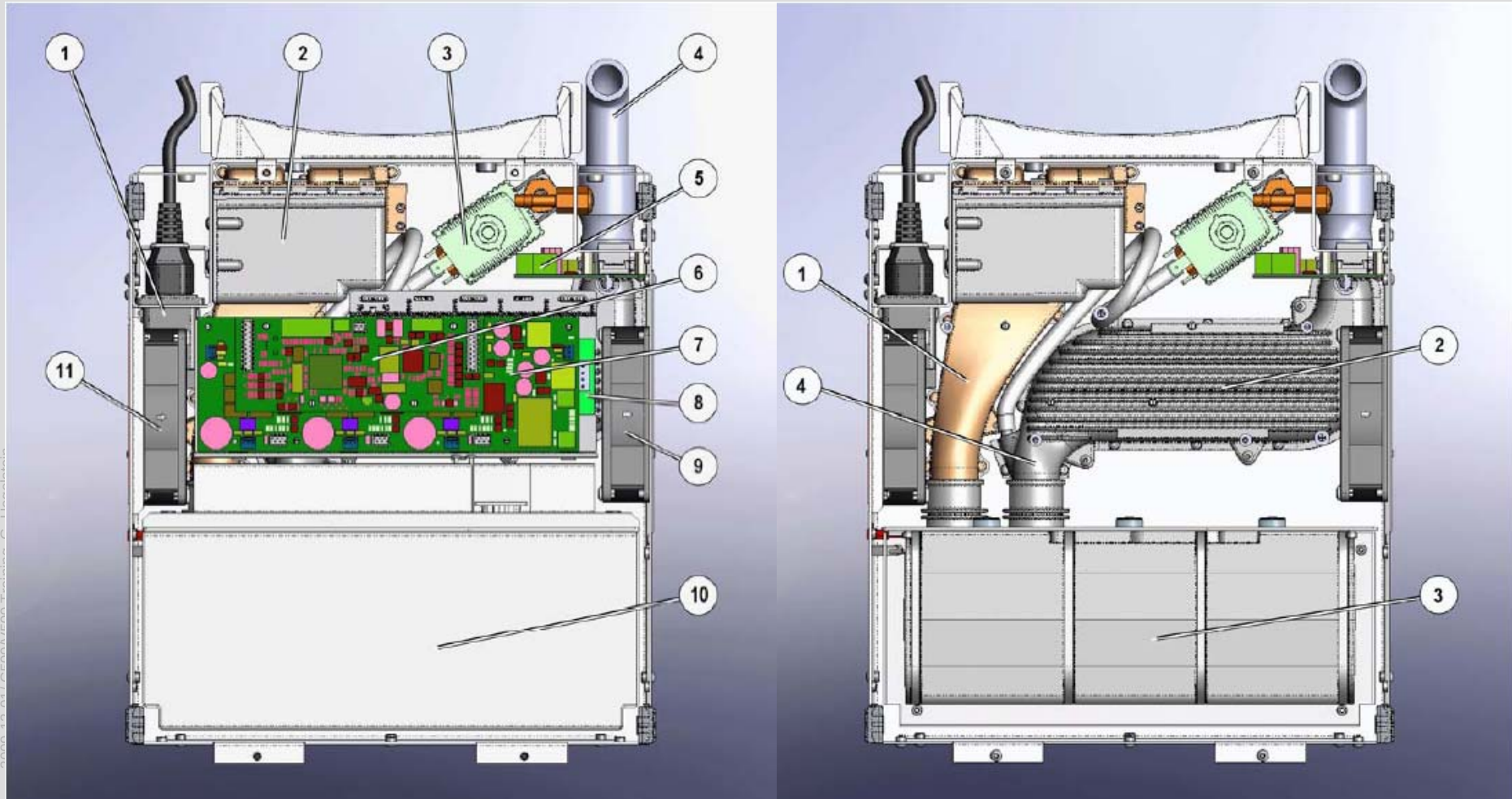

Workstation Critical Care / Neonatal Care / EvitaV300 Option GS500, function description

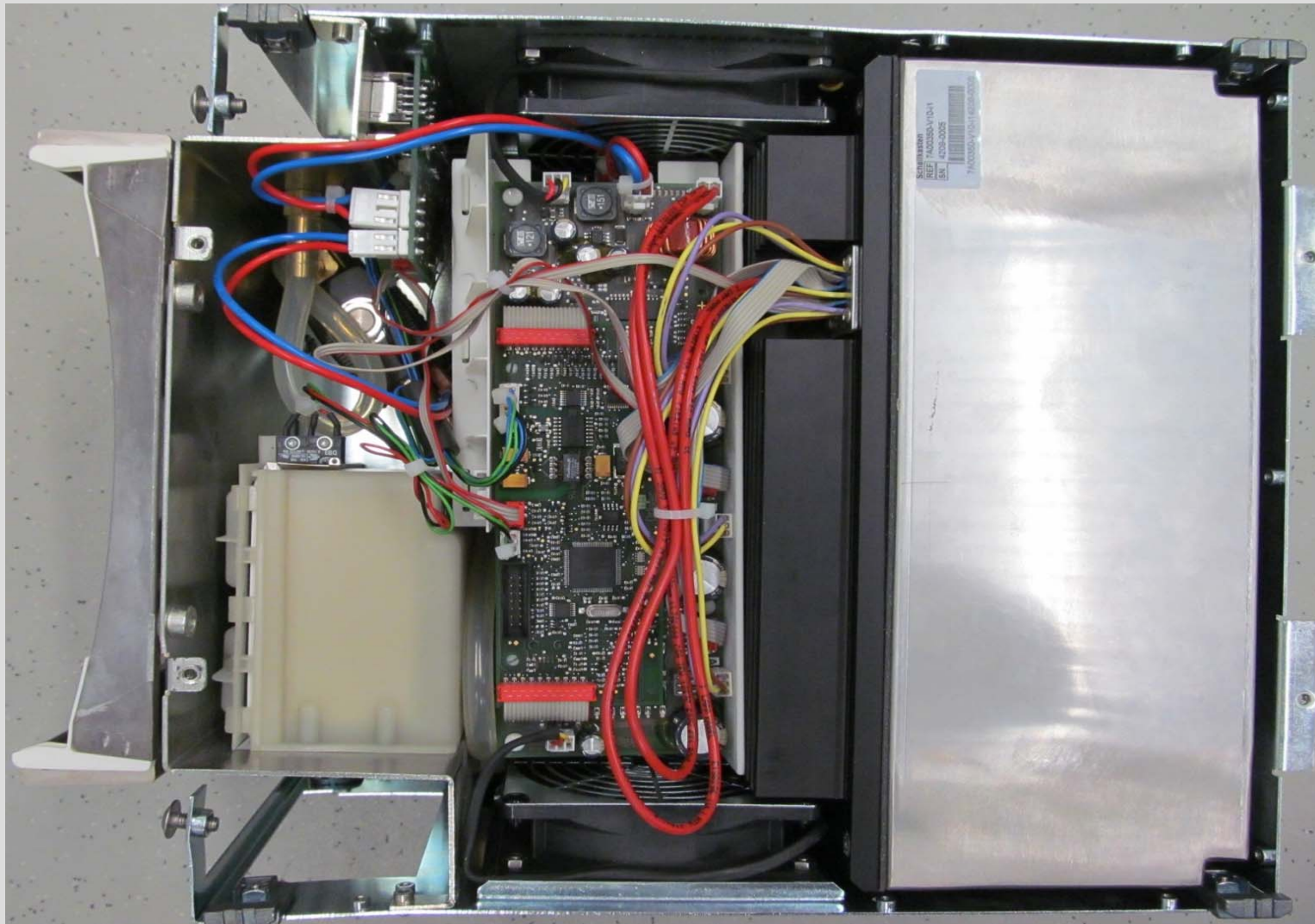


2009-12-01/ C500/V500 Training, G. Hagelstein

Workstation Critical Care / Neonatal Care / EvitaV300

Option GS500, function description



2009-12-01/ C500/V500 Training, G. Hagelstein

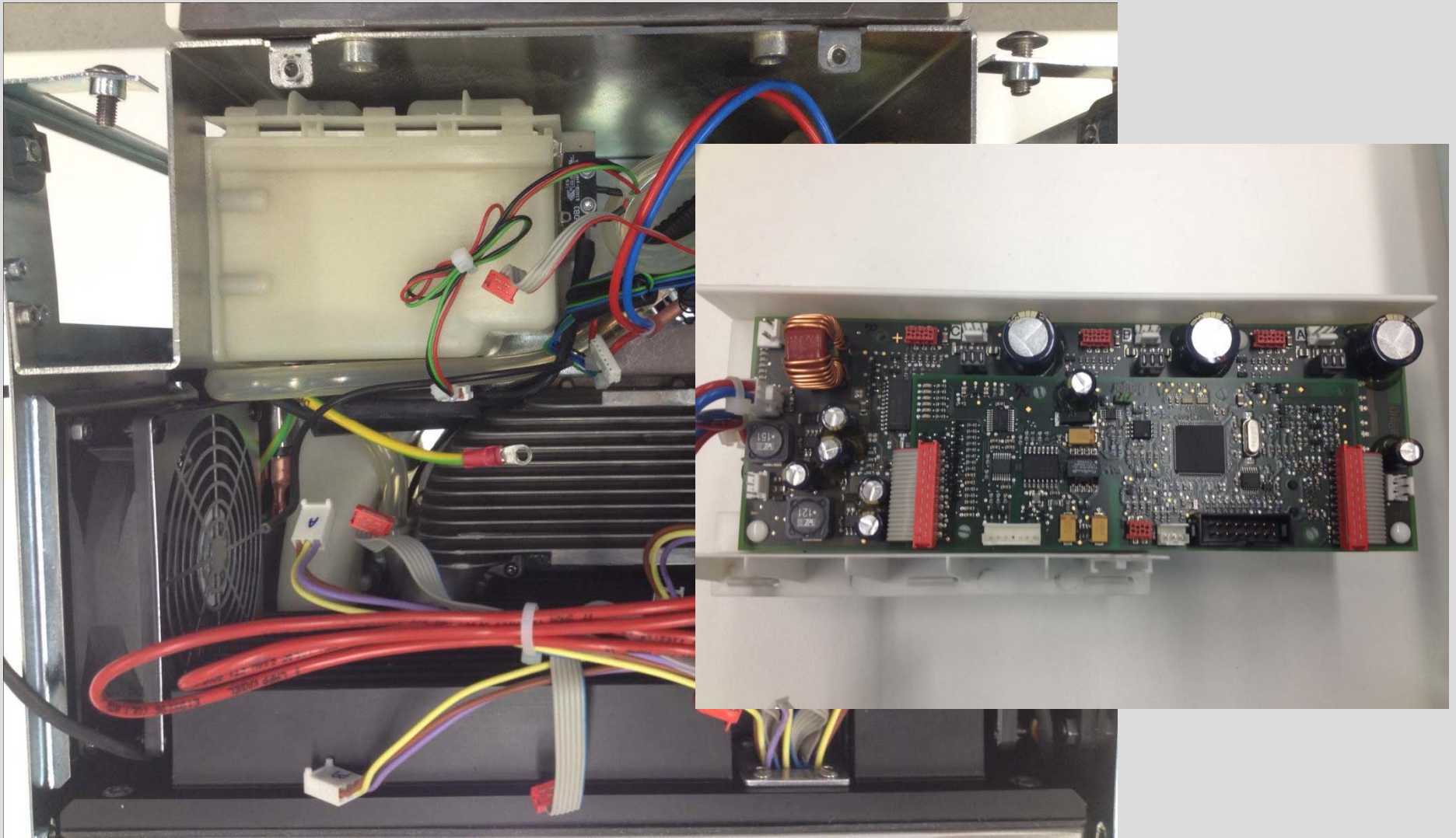

Workstation Critical Care / Neonatal Care / EvitaV300



Option GS500, function description



- Ambient air is drawn in via the respiratory gas filter from the blower unit and compressed.
- The speed of the radial fans in the blower unit is regulated by the controller.
- To cool the device and the respiratory gas, the ambient air is drawn into the housing by a fan.
- Sensors continuously measure the ambient temperature and humidity.
- The respiratory gas is cooled by way of a heat exchanger.
- The waste heat from the device and respiratory gas cooling is routed out of the housing by a fan.
- The respiratory gas is fed to the Evita Infinity V500 or Babylog VN500 ventilation unit.
- Data transfer and communication is routed via the serial port.

Workstation Critical Care / Neonatal Care / EvitaV300

Option GS500, technical data

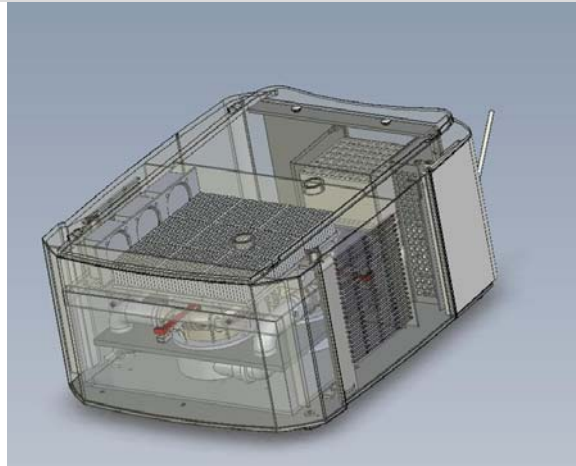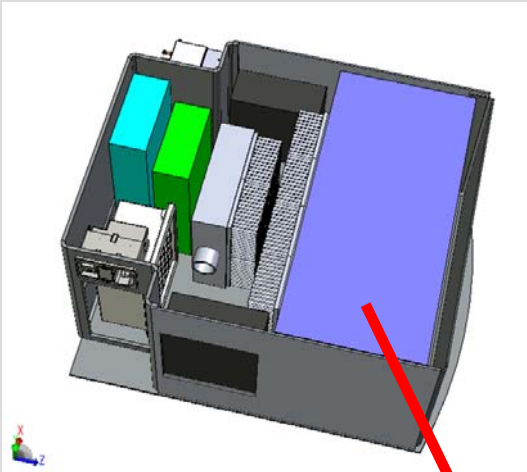


Technical Data:

ca. 120 mbar @ 100 lpm


ca. 80 mbar @ 180 lpm

Loudness = 48 dB(A)

Technology: 3-stage blower (Papst)

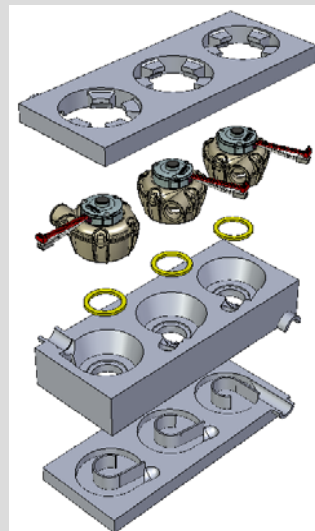

Workstation Critical Care / Neonatal Care / EvitaV300 Option GS500, function description blower unit

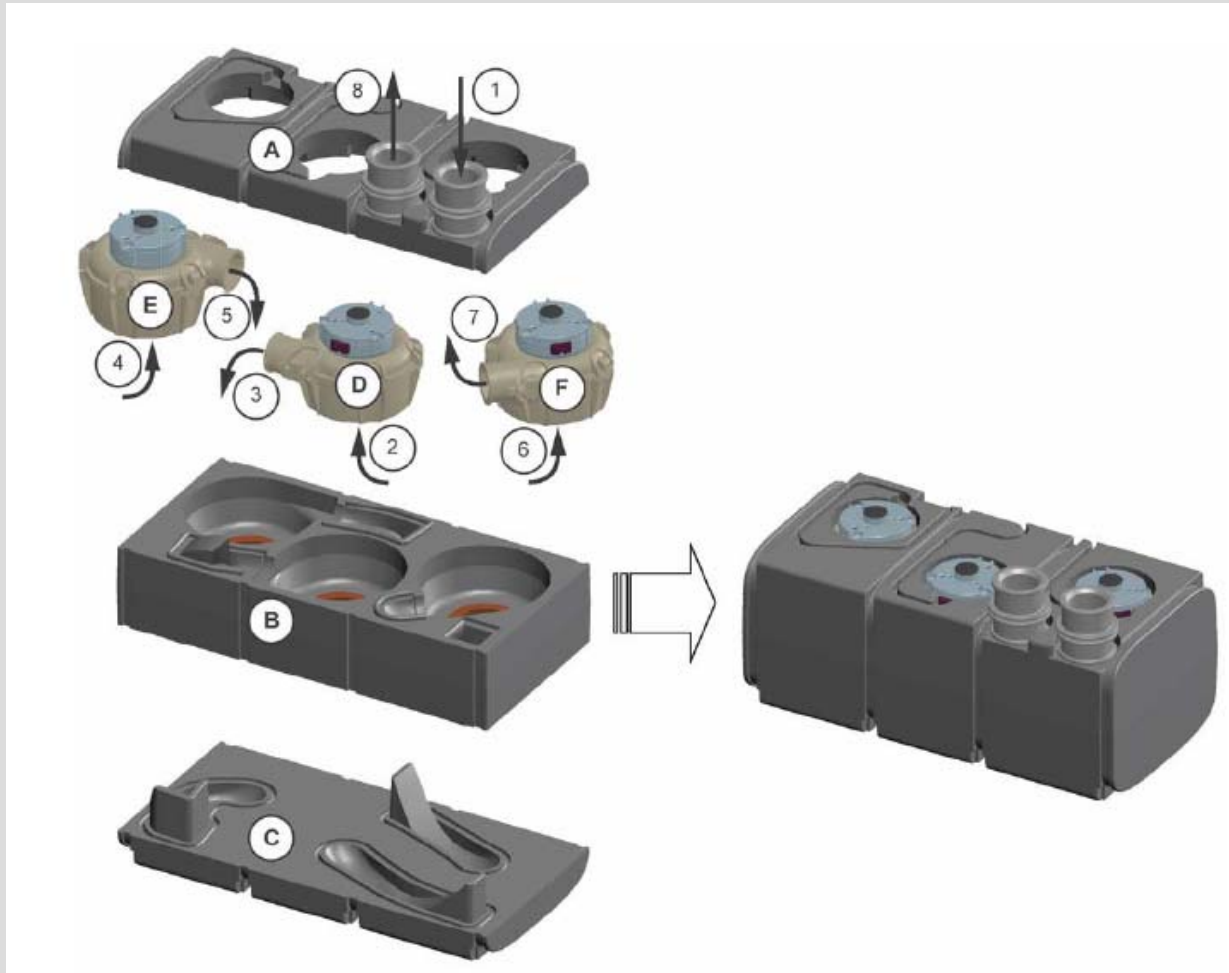

2009-12-01/ C500/V500 Training, G. Hagelstein



Workstation Critical Care / Neonatal Care / EvitaV300

Option GS500, software download



GS500

- **software download** via RS232 (Meditest) **direct on GS500**. Only possible with DrägerServiceSoftware 11.8 (or higher versions).
- device software for GS500 available in ServiceConnect; **actual is version 1.07** (published with TS215 SW2.50)

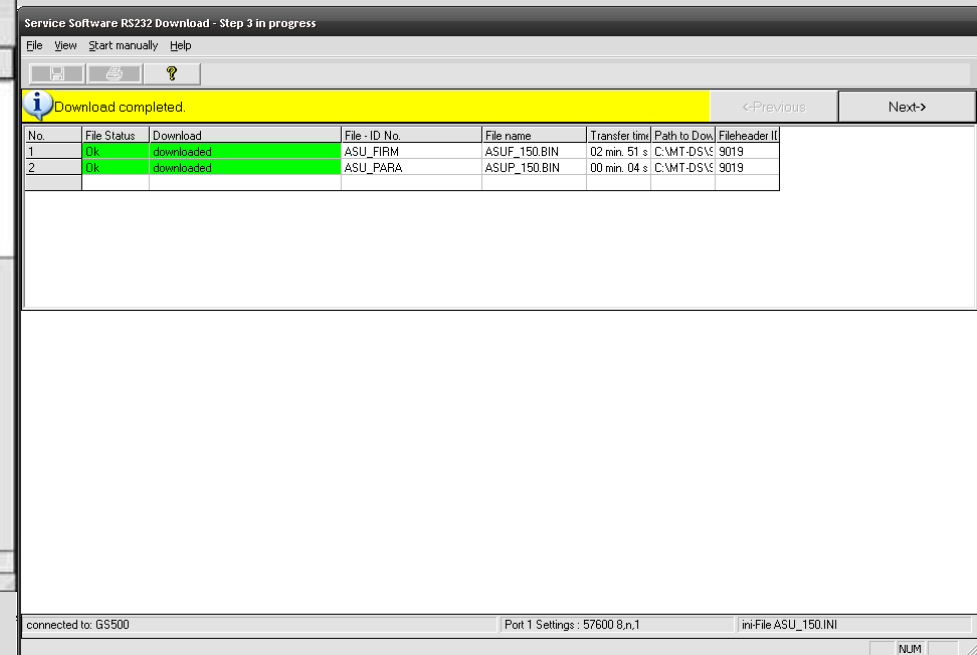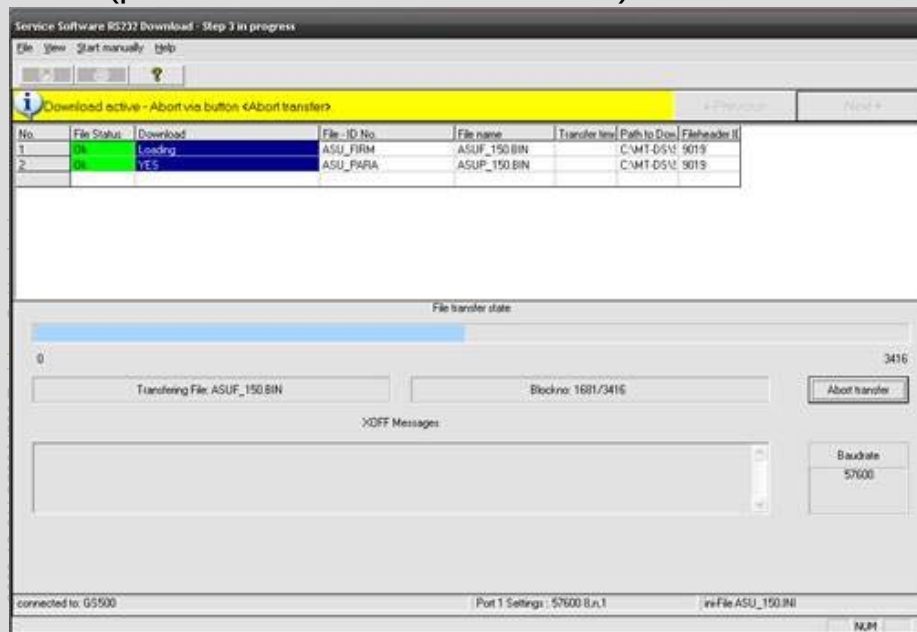

Workstation Critical Care / Neonatal Care / EvitaV300

Software GS500



Dräger

User: hagelstg Account number: DR HL Personal settings

Home Equipment Documents Parts **Software** Reports Remote Services Media

Product name
Evita V300

▼ Evita V300

- Evita V300 Ventilation unit
- GS500**
- Infinity C300 (MK31300)
- M11.1 MCable-Mainstream CO2 sensor

Software für GS500

VERSION FW 01.07
Softwareversion: GS500

This software update improves the system check related to the use of the GS500 at altitudes.
[...mehr](#)

[Herunterladen](#) [i](#)

VERSION FW 01.06.01
Softwareversion: GS500

Bugfix update:
[...mehr](#)

[Anfrage senden](#) [i](#)

2009-12-01/ C500/V500 Training, G. Hagelstein

GS500-DC

- GS500-DC is driven by an 24V-DC/DC-converter instead of an 220V power supply.



- Is the workplace running on main supply, will the GS500-DC get the power from power supply M7.3. When the main supply fails, the GS-500-DC get the power from the batteries of the workplace:
- Workplace with GS500-DC without PS500: the internal NiMH-battery of V-unit delivers 24V to GS500-DC. (Total running time of workstation is reduced from 30min to max. 15min!)
- Workplace with GS500-DC and PS500: the V-Unit doesn't have the internal battery. The batteries from PS500 are connected to M7.3. The GS500-DC is connected via a cable between M7.3 and PCB Power and will get the power via M7.3 from the batteries in the PS500. Running time: approx. 120min

GS500 – DC Installation

- Dräger has provided a conversion kit, with it an GS500-AC (230V) must be converted to GS500-DC (24V-DC/DC-converter).
- This conversion from GS500-AC (230V) to GS500-DC (24V) has to be done according TSB 08 (Class 2) until 31.12.2013!
- A workplace with the option „GS500-DC“ has a new PCB on M16 installed. This PCB (PBA Remote) is switching of the GS500-DC, when the workplace is shutting down. (The GS500-DC has no power switch.)

PCB PBA Remote

